

Accessories



Half Length Assist Device

Staff Control

Total Lock

Battery Back Up

Half Length Assist Device

Design is mounted at the perfect height to assist resident in and out of bed and provides continuous protection without the feeling of restraint for the resident.

Staff Control

Located within the footboard, the staff control gives caregivers easy access to bed operations. Features fluid protection rating and easy push buttons for easy operation and resident lockout.

Battery Back Up

Provides a low cost solution and protection for use during a power outage. Charger also available.

Total Lock

Provides four locking casters system for all 4 bed corners.



Wall Saver

Provides wall protection and reduces maintenance costs on headboards. Standard on most beds, ask for details.


Economy Lock Out (Not Shown)

An alternative to an embedded staff control, this lock out option locks out all functions on the pendant.



Assist Handle

Stable handhold assists with resident transfers in and out of bed and assists with resident repositioning on the sleep surface.

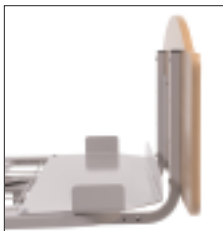

4" Bed Extension

The 4" extension is designed to extend the sleep surface an extra 4". This allows the use of the 80" mattress on a 76" deck and a 84" mattress on an 80" deck.



Side Saver

Side Saver provides protection down the length of the bed which is ideal for beds being repositioned frequently.




Trapeze and Trapeze Adapter

Trapeze adjusts to accommodate residents' reach. Trapeze adapter allows attachment of a trapeze system to a bed.



Wide Extension Kit

Transforms bed into a 42" sleep surface providing safer positioning, prevents resident falling out of bed and provides a more comfortable night's sleep.

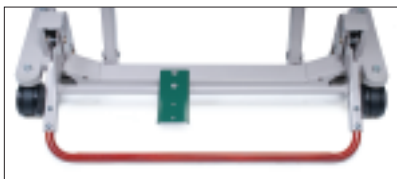

Z-Lok

Designed for mobility and security at any height with easy green "go" red "lock" system. Provides a one-step locking system that is easily accessible even with low bed height. Only available for the Zenith beds.



Pendant Holster

The pendant holder easily affixes to the assist device making access and storage simple. Easy read buttons make resident operation simple.



Edema Foot Ratchet

Provides quick positioning for improved circulation in the legs. Standard on most beds, ask for details.

Corner Mattress Retainer Kit

Retainers prevent mattress from sliding horizontally or vertically. Standard on most beds, ask for details.



EZ Lock/Flex Assist Devices

Rails provide resident protection to prevent rolling out of the bed. Rails slide down for easy entrance/exit from bed.

Panel IV Holder (Not Shown)

Holder secures on head/foot board eliminating the need for a separate standard IV stand.

For more information on Basic American accessories/options, contact your Basic American representative.

Note: Some accessories/options are not available for all bed models. Please contact your Basic American representative for details.



BASIC AMERICAN
MEDICAL PRODUCTS

800.554.9215

www.grahamfield.com