


★★★
PATRIOT
Homecare Bed Series

Introducing the new Patriot Full-Electric Homecare Bed with a GRID SLEEP SURFACE!

Designed and engineered by the same team who develop the GF Long Term Care beds.

Sleep surface is **MADE IN THE USA!** 

The single-motor design homecare bed is now available in a full-electric/low version. The Patriot Full-Electric Homecare Bed offers the same features as the Patriot Semi-Electric Homecare Bed with the additional feature of electric height adjustment. This makes it easier for caregivers to raise and lower the bed for patient care and linen changes. The electric height adjustment also enables the patient to easily position the bed for maximum comfort. The full-electric model can be adjusted to be a conventional height homecare bed OR a low bed for better inventory management!

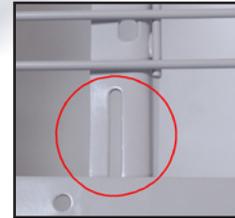

Slots in frame ensure proper rail placement



Locking pins secure frame to sleep deck

- Total assembled bed weighs only 178 pounds (bed and motor).
- Powerful DC motor for effortless height adjustment.
- Hi-Lo motor comes pre-installed on foot section for easy set-up and disassembly.
- All functions fully-integrated into convenient handheld pendant.
- Welded motor-guard protects motor and universal joint from damage during shipment and deliveries.
- Available with Fiberboard Decorative Walnut Bed-ends or High-impact Plastic Bed-ends.
- Designed to meet the requirements of applicable HCPCS codes.
- Maximum patient weight 350 lbs.
- 450 lb evenly distributed total weight capacity including patient, mattress, rails, and accessories.
- **Height Adjustment:**
Casters in Low Position:
 Min. 9.5", Max. 20" (at deck level)
Casters in Standard Position:
 Min. 13", Max. 23.5" (at deck level)

- Overall Bed Size: 87" L x 36" W.
- Lift Range: Head - 60°, Knee - 30°.
- No external pinch points!
- Slots in frame to ensure proper rail placement.
- Locking pins secure bed ends to sleep deck. Attached lanyards prevent loss while moving or storing.
- Only homecare bed to offer a 4-inch extension kit that can extend either the head, foot, or both head and foot decking. Bed length will be 84" with one kit, 88" with 2 kits.
- The only homecare bed to offer optional corner mattress retainers .

Grid Deck offers the following advantages:

- Much stronger frame.
- Better mattress support.
- No spring sag or bottoming out with heavier patients.
- Easier to clean.
- No connecting links or springs to lose.



Fiberboard Decorative Walnut Bed-Ends	High Impact Plastic Bed-Ends	Description	(UOM: 1 ea)	HCPCS
US0458	US0458PL	Full-Electric Bed Only		E0297
US0458-220V	US0458PL-220V	Full-Electric Bed Only - 220 V - International Only		
US0458-IFPKG	US0458PL-IFPKG	Full-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and Full Chrome Rails		E0265
US0458-IFPKGHRA	US0458PL-IFPKGHA	Full-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and Clamp-On Half Chrome Rails		E0265
US0458-PKG	US0458PL-PKG	Full-Electric Bed with 1633-Innerspring Mattress and Full Chrome Rails		E0265
US0458-PKGHRA	US0458PL-PKGHRA	Full-Electric Bed with 1633-Innerspring Mattress and Clamp-On Half Chrome Rails		E0265
US0458-RFPKG	US0458PL-RFPKG	Full-Electric Bed with 1633-Rolled Foam Mattress and Full Chrome Rails		E0265
US0458-RFPKGHRA	US0458PL-RFPKGHA	Full-Electric Bed with 1633-Rolled Foam Mattress and Clamp-On Half Chrome Rails		E0265
US0458-RPKG	US0458PL-RPKG	Full-Electric Bed with No Mattress and Full Chrome Rails		E0265
US0458-RPKGHRA	US0458PL-RPKGHRA	Full-Electric Bed with No Mattress and Clamp-On Half Chrome Rails		E0265
US0458-XFIPKGF	US0458PL-XFIPKGF	Full-Electric Bed with 1633-Extra Firm Innerspring Mattress and Full Chrome Rails		E0265
US0458-XFIPKGHRA		Full-Electric Bed with 1633-Extra Firm Innerspring Mattress and Clamp-On Half Chrome Rails		E0265

Warranty:

- Two year warranty for defects in workmanship and materials of mechanical components, frame, and electronics.
- During the warranty period, defective items will be repaired or replaced at manufacturer’s option. The warranty does not include any labor charges incurred in replacement part(s) installation or any associated freight or shipping charges to the manufacturer.

Bed shipped in three cartons:

1. Foot section and motor.
2. Head section and Hi/Lo rod.
3. Bed end set.



Fiberboard Decorative Walnut Bed-Ends	High Impact Plastic Bed-Ends	Description	(UOM: 1 ea)	HCPCS
US0458	US0458PL	Full-Electric Bed Only		E0297
US0458-220V	US0458PL-220V	Full-Electric Bed Only - 220 V - International Only		
US0468-IFPKGFR	US0468PL-IFPKGFR	Full-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and FDA Full Rails		E0265
US0468-IFPKGHR	US0468PL-IFPKGHR	Full-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and FDA Half Rails		E0265
US0468-IFPKGQR	US0468PL-IFPKGQR	Full-Electric Bed with 1633-Innerspring Mattress with Foam Sleep Surface and FDA Quarter Rails		E0265
US0468-PKG	US0468PL-PKG	Full-Electric Bed with 1633-Innerspring Mattress and FDA Full Rails		E0265
US0468-PKGHR	US0468PL-PKGHR	Full-Electric Bed with 1633-Innerspring Mattress and FDA Half Rails		E0265
US0468-PKGNR	US0468PL-PKGNR	Full-Electric Bed with 1633-Innerspring Mattress and No Rails		E0265
US0468-PKGQR	US0468PL-PKGQR	Full-Electric Bed with 1633-Innerspring Mattress and FDA Quarter Rails		E0265
US0468-RFPKG	US0468PL-RFPKG	Full-Electric Bed with 1633- Rolled Foam Mattress and FDA Full Rails		E0265
US0468-RFPKGHR	US0468PL-RFPKGHR	Full-Electric Bed with 1633- Rolled Foam Mattress and FDA Half Rails		E0265
US0468-RFPKGQR	US0468PL-RFPKGQR	Full-Electric Bed with 1633- Rolled Foam Mattress and FDA Quarter Rails		E0265
US0468-RPKG	US0468PL-RPKG	Full-Electric Bed with No Mattress and FDA Full Rails		E0265
US0468-RPKGHR	US0468PL-RPKGHR	Full-Electric Bed with No Mattress and FDA Half Rails		E0265
US0468-RPKGQR	US0468PL-RPKGQR	Full-Electric Bed with No Mattress and FDA Quarter Rails		E0265
US0468-XFIPKGFR	US0468PL-XFIPKGFR	Full-Electric Bed with Extra Firm Innerspring Mattress and FDA Full Rails		E0265
US0468-XFIPKGHR	US0468PL-XFIPKGHR	Full-Electric Bed with Extra Firm Innerspring Mattress and FDA Half Rails		E0265
US0468-XFIPKGQR	US0468PL-XFIPKGQR	Full-Electric Bed with Extra Firm Innerspring Mattress and FDA Quarter Rails		E0265

Warranty:

- Two year warranty for defects in workmanship and materials of mechanical components, frame, and electronics.
- During the warranty period, defective items will be repaired or replaced at manufacturer’s option. The warranty does not include any labor charges incurred in replacement part(s) installation or any associated freight or shipping charges to the manufacturer.

Bed shipped in three cartons:

1. Foot section and motor.
2. Head section and Hi/Lo rod.
3. Bed end set.