



**LUMEX** ★★★  
**PATRIOT<sup>TX</sup>**  
 HOMECARE BED

Graham-Field carries a complete line of mattresses and bed accessories to provide the best quality and comfort for your homecare and hospice needs.

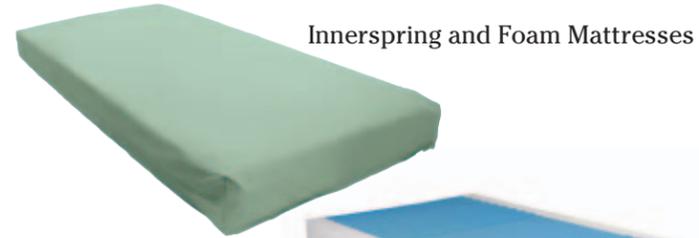

Innerspring and Foam Mattresses

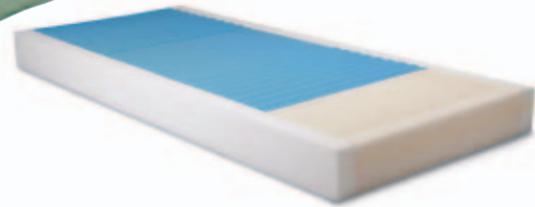

Trapezes



Overbed Tables



Hydraulic and Power Patient Lifts



Patient Alarms and Floor Mats



Bed Accessories

And More!



Manufactured by:  
 GF Health Products, Inc.  
 Fond du Lac, WI 54937  
*Proudly made in the USA*  
 770.368.4700 [www.grahamfield.com](http://www.grahamfield.com)

© 2014, GF Health Products, Inc. GF1300084RevB14

**LUMEX** ★★★  
**PATRIOT<sup>TX</sup>**

HOMECARE BED  
*PROUDLY MADE IN THE U.S.A.*



The Lumex Patriot LX Homecare bed is designed to be the most advanced and quality driven homecare bed in the market. Weighing only 118 lbs and with a unique and attractive design to blend in the home, the Patriot LX now offers a lighter weight and all encompassing solution for your homecare and hospice needs.

**LUMEX** ★★★  
**PATRIOT LX**  
 HOMECARE BED



**Lightweight  
118 lbs!**




Robust grid deck sleep surface provides superior mattress support and easy to clean surface compared to traditional spring decks. It features **SilverSolutions**- a anti-microbial technology that protects against bacteria and odor.



Patent pending integrated hi/lo mechanism allows for smooth and easy bed height adjustment and reduces potential loss of parts such as high-low rod when used as rentals.



Self contained crank for semi-electric can be used at either end to make adjustment from low 10" to 20" easy and even allows for Trendelenburg and slant positioning. The full-electric model features dual hi-lo motors to adjust bed height.

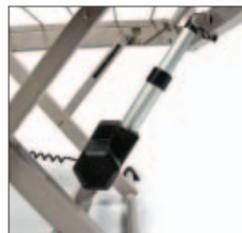

Single motor unit for head and foot articulation is designed for easy installation, requires less maintenance and provides the ultimate in quiet, smooth operation.



Attractive, "home-like" thermofused bed ends offer a less institutional look and eliminate a bed component that often requires frequent repair and maintenance.

Other Features include:

- Lightweight, patent pending design reduces the risk of injuries for technician during setup and storage.
- Head and foot sections can be raised or lowered simultaneously or independently with one hand control.
- During a power failure, the motor can lower both the head and foot section with just a single 9 volt battery.
- 4" extension kit can be installed on head and/or foot to extend bed length to 84" or 88".

US5000 Semi-Electric

**The Cost of Quality**

- Reduced cost - the Patriot LX is built to reduce the overall cost of ownership because it is designed to last over time
- Reduced risks - weighing only 118 lbs (126 lbs for full-electric) reduces the amount of back injuries to your bed technicians and property damage during transport and setup
- Increased efficiency - bed setup is so simple that it can save on time and requires minimal maintenance keeping your drivers time efficient



**Liberty Rails**

Lumex Patriot™ Liberty Bed Rails have been designed to offer a safe and secure bed rail system to provide you with the utmost level of quality to meet your homecare bed needs. Developed using the **FDA Hospital Bed System Dimensional Guidance to Reduce Entrapment** and tested to the **Hospital Bed Safety Workgroup and FDA Dimensional Test Methods for Bed System Zones 1-4 and International Standard IEC60601-2-38 Requirements for the Safety of Electrically Operated Hospital Beds**, these new FDA Compliant bed rail systems offer added protection and reduce the risk of entrapment.



**Lumex® Full Length Rail (GF6570B-1)**

The Full Length Rail is a nominal, telescoping full length rail that provides maximum length rail perimeter coverage, intended for use by patients requiring maximum length rail and maximum perimeter protection.

**Lumex® Half No-Gap Head Rail (GF6590BH-1)**

Provides full head coverage of 31.5 inches and extends to eliminate the gap between the head board and the end of the rail.




**Features include:**

- Improved rail design decreases the spacing between the crossbars, thus reducing the risk of entrapment. New bed rail length closes open areas that can occur during full bed articulation.
- Rails now mount in fixed patent pending rail slots (shown below) on Patriot™ bed deck to prevent sliding and moving — easy installation saves delivery and setup time and ensures proper rail placement.
- Spring-loaded, self-adjusting crossbraces feature safety tabs that prevent accidental disengagement when bumped or during application of side pressure.



Patent pending rail slot

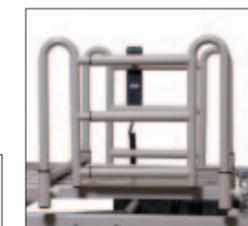

**Lumex® Quarter Rail (GF6580BH-1)**

Approximately 18 inches in length when mounted, this rail can be used as a head assist, head rail, or foot rail to provide minimal rail coverage and gives assistance when entering or exiting the bed.

**Technical Specifications**

- Maximum patient weight 350 lbs **EVENLY DISTRIBUTED**
- 450 lb distributed total weight capacity including patient, mattress, rails and accessories
- Overall Bed Size: 84" x 35"
- Lift Range: Head – 60°, Knee – 30°
- Bed Height Range: 10" - 20" - no conversion kit needed!
- Product weight: US5000 Semi-Electric 118 lbs, US6000 Full-Electric 126 lbs

Warranty: Frame – 5 Years, Electronics – 2 Years, All Other Components – 1 Year