

The Traveler® HD is manufactured with the same commitment to quality E&J is known for.

Designed for optimal patient transport as well as the hospital/clinical use. Built with a heavy duty frame to support up to 500 lb users. Reinforced at all critical load-bearing points and a long lasting embossed back and seat upholstery with a double inner liner for durability. With its multiple seat widths ranging from 20" to 24" and 18" seat depth (1" longer than standard), this heavy duty wheelchair enables institutions and caregivers to accommodate a large variety of patients.

### Standard Features





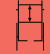


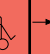
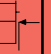

- Available in 20", 22" and 24" seat widths, with 18" seat depth (1" longer than standard)
- Carbon steel reinforced frame include double cross braces, heavy-duty 5/8" diameter axles
- Embossed black leatherette upholstery with double inner lining for durability with standard chart pocket on back
- Custom color program is available in Grey, Orange or Burgundy vinyl upholstery for identification control
- Hammertone paint provides a long lasting and durable finish
- Two-position axle and caster fork provide 2" seat height adjustment for a hemi position
- Tool-free adjustable front riggings with durable aluminum black footplates
- Rear 24" x 1" black plastic wheels with black handrims and non-marking polyurethane tires
- Maximum Weight Capacity 500 lb - EVENLY DISTRIBUTED
- Meets Medicare code K0007



### Warranty

Five (5) year limited warranty



<b>Traveler® HD</b>		<i>Dimensions and Weights</i>									
		<i>Seat</i>			<i>Back</i>	<i>Overall</i>				<i>Weight</i>	
		<i>Width</i>	<i>Depth</i>	<i>Height</i>	<i>Height</i>	<i>Depth</i>	<i>Height</i>	<i>Width</i>	<i>Folded Width</i>	<i>Without front rigging</i>	<i>Max Capacity</i>
											
<b>3G010320</b>	Detachable Desk Arm, Swingaway Footrest	20"	18"	18" or 20"	16 ½"	34"	36"	28 ½"	12 ½"	57 lb	500 lb
<b>3G010330</b>	Detachable Desk Arm, Elevating Legrest	20"	18"	18" or 20"	16 ½"	34"	36"	28 ½"	12 ½"	57 lb	500 lb
<b>3G010340</b>	Detachable Full Arm, Swingaway Footrest	20"	18"	18" or 20"	16 ½"	34"	36"	28 ½"	12 ½"	57 lb	500 lb
<b>3G010350</b>	Detachable Full Arm, Elevating Legrest	20"	18"	18" or 20"	16 ½"	34"	36"	28 ½"	12 ½"	57 lb	500 lb
<b>3G010420</b>	Detachable Desk Arm, Swingaway Footrest	22"	18"	18" or 20"	16 ½"	34"	36"	30 ½"	12 ½"	57 lb	500 lb
<b>3G010430</b>	Detachable Desk Arm, Elevating Legrest	22"	18"	18" or 20"	16 ½"	34"	36"	30 ½"	12 ½"	57 lb	500 lb
<b>3G010440</b>	Detachable Full Arm, Swingaway Footrest	22"	18"	18" or 20"	16 ½"	34"	36"	30 ½"	12 ½"	57 lb	500 lb
<b>3G010450</b>	Detachable Full Arm, Elevating Legrest	22"	18"	18" or 20"	16 ½"	34"	36"	30 ½"	12 ½"	57 lb	500 lb
<b>3G010520</b>	Detachable Desk Arm, Swingaway Footrest	24"	18"	18" or 20"	16 ½"	34"	36"	32 ½"	12 ½"	57 lb	500 lb
<b>3G010530</b>	Detachable Desk Arm, Elevating Legrest	24"	18"	18" or 20"	16 ½"	34"	36"	32 ½"	12 ½"	57 lb	500 lb
<b>3G010540</b>	Detachable Full Arm, Swingaway Footrest	24"	18"	18" or 20"	16 ½"	34"	36"	32 ½"	12 ½"	57 lb	500 lb
<b>3G010550</b>	Detachable Full Arm, Elevating Legrest	24"	18"	18" or 20"	16 ½"	34"	36"	32 ½"	12 ½"	57 lb	500 lb
<b>3G020500</b>	Fixed Full Arm, Fixed Footrest	24"	18"	18" or 20"	16 ½"	34"	36"	32 ½"	12 ½"	57 lb	500 lb

<b>Traveler® HD – Common Accessories &amp; Parts</b>		<b>HCPSC Code</b>
<b>90763064</b>	Heel loop, pair	E0951
<b>90763060</b>	Wheel lock, Push-to-lock (includes rubber tips) pair	E2206
<b>90763061</b>	Wheel lock, Pull-to-lock (includes rubber tips) pair	E2206
<b>90763062</b>	Wheel lock extension, pair	E0961
<b>90763063</b>	Rubber tip, wheel lock (10 each)	
<b>300-505</b>	Positioning belt (auto-style) push button 27" to 47"	E0978
<b>92133040</b>	Positioning belt (airline style) (metal buckle) 36"	E0978
<b>90133125</b>	Positioning belt long (airline style) (metal buckle) 47"	E0978
<b>B10321G</b>	O2-IV holder for HTC	E2208
<b>B10325</b>	O2 cylinder holder	E2208
<b>90763065</b>	Anti-tippers, pair (includes hardware) for 20" seat height Traveler SE & HD	E0971
<b>90763066</b>	Anti-tippers, pair (includes hardware) for 18" seat height Traveler SE & HD	E0971
<b>92523850</b>	Anti-fold bar (chromed) lower frame for Advantage, SE, LX, L3, L4, HD & HTC	