



CLINICAL CARE RECLINERS

PIVOT-ARM CLINICAL CARE RECLINER

Three reclining positions for patient comfort. Plus a caregiver-activated Trendelenburg (or Shock) Position, is easily actuated from either side by a conveniently located foot pedal.

Both arms pivot outward, enabling both cleaning and easier patient transfers.

Wide seat ensures comfort for different patient types and sizes.

Two self storing Side Tables for patient prep/convenience; smooth surface allows for easy cleaning.

Four Tente[®] swivel casters facilitate transport and maneuverability in confined areas. Two fixed front locking casters, one directional (swivel) locking rear caster, and one total locking rear caster.

Urethane armrest caps are comfortable for patient and easy to maintain; also available with fully upholstered arms.

Optional extended legrest provides comfort and better leg support for taller persons. *Extended legrest models do not include pull out foot plate*


All Recliner components can be replaced greatly enhancing the lifespan of the recliner and reducing its cost in use over time.

Wide selection of fabrics and colors. See the [Swatch Card](#) for details.



Options:

- Heat and Massage
- Heat Only feature
- Fully Upholstered Armrests
- Extended Legrest
- IV Pole Holder / IV Pole
- Headrest Cover
- Universal Pillow

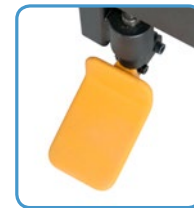

Multiple Positions for Enhanced Patient Comfort



GF Health Products, Inc. is an
intertek ISO 13485:2016 / MDSAP Certified Company

Color variations may occur due to monitor or printer settings.

+1 770.368.4700

Information contained herein is subject to change without notice.
The most current and complete product information can be found on our website.

www.grahamfield.com



© 2019, GF Health Products, Inc. All Rights Reserved.
Lumex is a registered trademark of GF Health Products, Inc.



CLINICAL CARE RECLINERS

PIVOT-ARM CLINICAL CARE RECLINER



Model (xxxx represents fabric colors)	Wall Saver Armrest Edges	Fully Upholstered Arms	Heat & Massage	Heat Only	Extended Legrest	overall Width		overall Depth		overall Height		seat Width		seat Depth		seat Height		Arm Height from Seat		Width between arms		Width between side tables in upright position		Depth fully reclined (Trendelenburg)		Weight Capacity EVENLY DISTRIBUTED		UOM
						inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	lb	Kg	
FR601P	●					34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	72	182.9	400	181.4	1 ea
FR601PH	●	●				34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	72	182.9	400	181.4	1 ea
FR601PHO	●			●		34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	72	182.9	400	181.4	1 ea
FR601UxxxxPA		●				34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	72	182.9	400	181.4	1 ea
FR601UHxxxxPA		●	●			34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	72	182.9	400	181.4	1 ea
FR601UHOxxxxPA		●		●		34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	72	182.9	400	181.4	1 ea
FR601ULxxxx	●				●	34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	76	193.1	400	181.4	1 ea
FR601UHLxxxx	●	●			●	34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	76	193.1	400	181.4	1 ea
FR601UHOxxxx	●			●	●	34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	76	193.1	400	181.4	1 ea
FR601PLxxxx		●			●	34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	76	193.1	400	181.4	1 ea
FR601PHLxxxx		●	●		●	34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	76	193.1	400	181.4	1 ea
FR601PHOLxxxx		●		●	●	34.5	87.6	38	96.5	51.5	130.8	22.5	57.2	19	48.3	22	55.9	30	76.2	25	63.5	54	137.2	76	193.1	400	181.4	1 ea



GF Health Products, Inc. is an
intertek ISO 13485:2016 / MDSAP Certified Company

Color variations may occur due to monitor or printer settings.

+1 770.368.4700

Information contained herein is subject to change without notice.
The most current and complete product information can be found on our website.

www.grahamfield.com



© 2019, GF Health Products, Inc. All Rights Reserved.
Lumex is a registered trademark of GF Health Products, Inc.