



EXTENDED CARE BEDS



MATRIX SERIES

The Matrix bed features sealed ball bearings in the joints for smooth, quiet movement and SilverSolutions in the powder coating to provide antimicrobial surface protection.

Basic American by Graham-Field is a leading manufacturer of beds serving the senior living industry for 40+ years. We proudly manufacture our beds in the USA, in Fond du Lac, Wisconsin, which enables faster response to customers' needs.

Innovation with industry-leading technology is designed into the Basic American Zenith, Matrix, and Liberty models to provide feature-rich beds and value for every budget.



6100
MATRIX BED
*shown in chair position with
optional two-pedal locking system*

MATRIX SERIES

FEATURES AND ADVANTAGES:

Made in the USA

- Minimal lead times
- Quality craftsmanship and manufacturing

SilverSolutions

- Located in the powder-coating which provides antimicrobial surface protection
- Premium Products

Sealed Ball Bearings

- Allows seamless, continuous, and noiseless movement
- Maintenance Free

Expandable Deck Available

- Pin-on wide kits extend the deck width to 39" or 42"

ANSI/AAMI ES60601 Safety Certification

- Tip stability testing
- Load capacity testing (twice the operating capacity)
- Roll clearance testing
- Splash proof electronics testing

Optional 'Plug and Play' digital electronics enable future feature additions

- Staff Control for Trendelenburg and Reverse Trendelenburg
- Chair Position
- Underbed Light
- USB Charger

Certified as a transport device for the CMS Emergency Preparedness Plan

- Easy maneuvering/adjustability through doorways and hallways during emergency or evacuation

Optional Intuitive two-pedal locking system

- Centrally located
- Easily visible and accessible at ALL heights and widths

Industry-leading, energy-saving electronics

- Consume up to 66% less power than transformer electronics

Grounded bed with right-angle hospital grade plug

- Provides added safety for both staff and residents

US Patented design of high/low mechanism

- Housed with the mainframe tubes; no holes are located on the top side that will harbor bacteria


Grid sleep surface with cross braces

- Sturdy, no-sag sleeping surface



5100
MATRIX BED

USB charger (right)
underbed light (below)

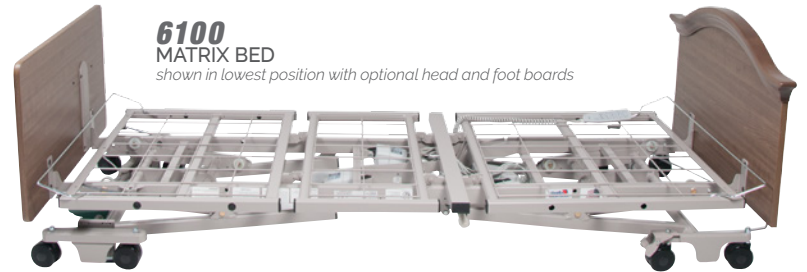

6100
MATRIX BED

shown in lowest position with optional head and foot boards

Warranty:

180	Main Frame & Welds ¹
48	Control Box & Actuator Motors
36	Hand Control Pendant
36	Staff Control & Cabling
12	Head Board & Foot Board
24	All Other Durable Components not otherwise listed


¹ 15 year frame warranty



6100

MATRIX BED

shown in the Trendelenburg position and with optional pedal lock

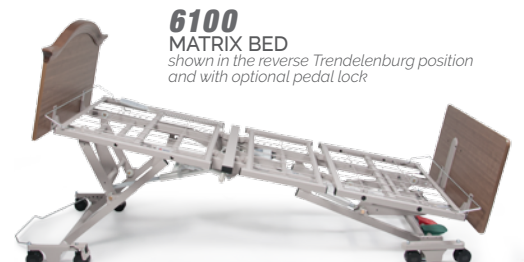

6100

MATRIX BED

shown in the reverse Trendelenburg position and with optional pedal lock



GF Health Products, Inc. is an
intertek ISO 13485:2016 / MDSAP Certified Company

Color variations may occur due to monitor or printer settings.

+1 770.368.4700

Information contained herein is subject to change without notice.
The most current and complete product information can be found on our website.
www.grahamfield.com



© 2021, GF Health Products, Inc. All Rights Reserved.
Basic American Medical Products is a trademark of GF Health Products, Inc.

graham
field