



## SUPPORT SURFACES





# GOLD CARE FOAM MATTRESS

## 416 & 419 SERIES


The Lumex 419 Series Foam Mattress provides support to the patient for pressure redistribution. The Lumex mattresses are specifically developed with patient comfort and pressure injury management in mind. Provides pressure reduction and redistribution.

Position	average Pressure mmHg
Scapula	14
Sacral Prominence	20
Heel	23
Trochanter	37

Tissue Interface Pressure Map: Gold Care Plus<sup>2</sup>

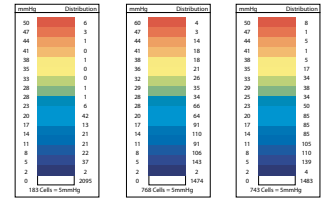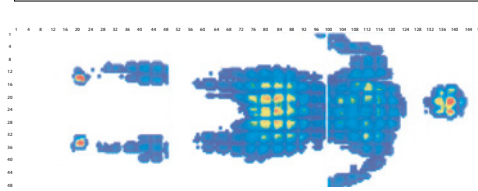

data on file :: 419 Gold Mattress  
June 10, 1999 on Male - 5' 6" - 165 lb

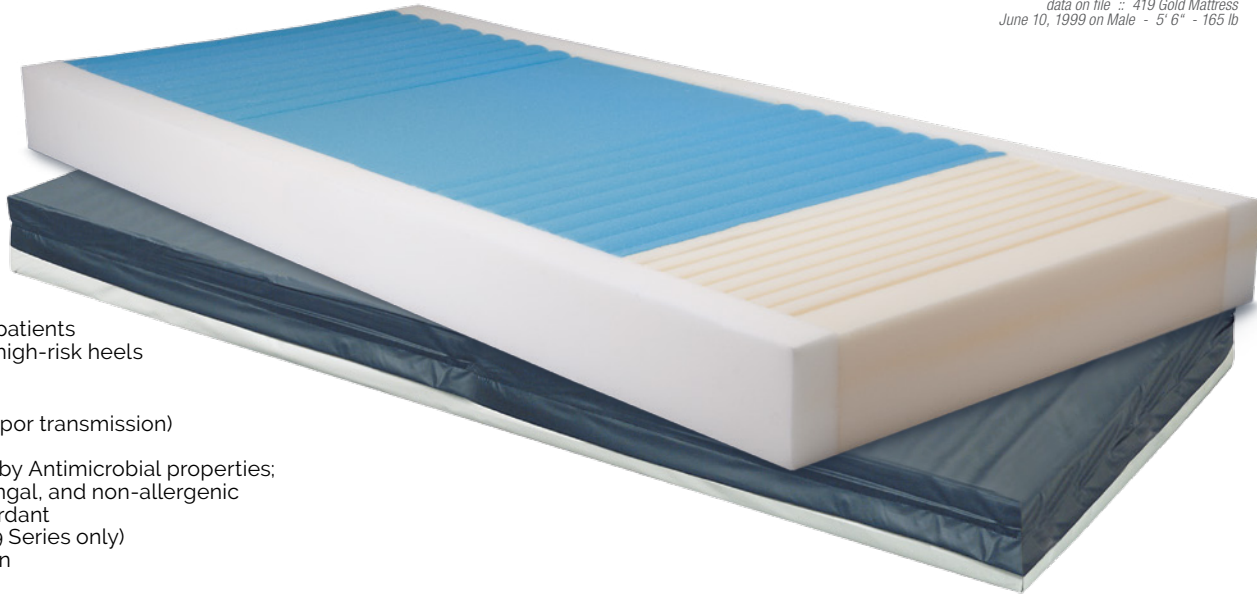

### Features

- Protection for higher risk patients with pressure ulcers and high-risk heels

### Cover features

- Vapour permeable (0% vapor transmission)
- Fluid resistant
- Infection Control is aided by Antimicrobial properties; Bacteria resistant, anti-fungal, and non-allergenic
- Anti-static and flame retardant
- Concealed full zipper (419 Series only)
- Easy to clean and maintain

### Mattress Construction

- Four separate anatomically matched support zones enhance patient comfort, reduce tissue interface pressures, and providing pressure redistribution
- Die-cut design on foam top provides comfort and allows air to flow and cool the sleep surface
- Firm side perimeter provides solid edge support and facilitates patient transfers
- Designed to meet the requirements of HCPCS Code E0184

### Specifications

Item	UOM	Meets HCPCS	Overall Width		Overall Length		Overall Thickness		Product Weight <sup>1</sup>		Weight Capacity EVENLY DISTRIBUTED		Raised Side Bolsters	Defined Perimeter	Soft Heel Zone	Soft Heel Slope	Pressure Reduction	Zippered Cover	Meets 16CFR1633	Made in USA	Mattress in months	Cover in months
			inch	cm	inch	cm	inch	cm	lb	kg	lb	kg										
41680-1633	1/ea	E0184	35	88.90	80	203.20	6	15.24	20.5	9.30	400	181.44	—	■	■	—	■	—	■	■	60	12
41980-1633	1/ea	E0184	35	88.90	80	203.20	6	15.24	20.5	9.30	400	181.44	■	■	■	■	■	■	■	■	60	12
41984-1633	1/ea	E0184	35	88.90	84	213.36	6	15.24	20.5	9.30	400	181.44	■	■	■	■	■	■	■	■	60	12
4168042-1633	1/ea	E0184	42	106.68	80	203.20	6	15.24	20.5	9.30	400	181.44	■	■	■	■	—	■	■	■	60	12
4198042-1633	1/ea	E0184	42	106.68	80	203.20	6	15.24	20.5	9.30	400	181.44	■	■	■	■	■	■	■	■	60	12
4198442-1633	1/ea	E0184	42	106.68	84	213.36	6	15.24	20.5	9.30	400	181.44	■	■	■	■	■	■	■	■	60	12

■ Standard ■ Optional ■ Available ■ Not Available ■ Accessory — Does Not Apply ● Active Hyperlink  
<sup>1</sup> approximate - as size changes so does weight to a minor degree  
<sup>2</sup> data on file



GF Health Products, Inc. is an  
 Intertek ISO 13485:2016 / MDSAP Certified Company

Color variations may occur due to monitor or printer settings.

Product Hyperlinks Enabled  
 +1 770.368.4700

Information contained herein is subject to change without notice.  
 The most current and complete product information can be found on our website.

[www.grahamfield.com](http://www.grahamfield.com)



© 2007, GF Health Products, Inc. All Rights Reserved.  
 Lumex is a registered trademark of GF Health Products, Inc.

graham  
 field

GF0700110 Rev H22