



# EXTENDED CARE BEDS



We proudly make all of our beds in the USA to the highest safety standards:<sup>4</sup>

- FDA Recognized Standard: ANSI/AAMI STD ES60601-1
- Health Canada Recognized Standard: CAN/CSA C22.2 No. 60601-1 (IEC 60601-1:2012-Edition 3.1)
- Includes International Standards: IEC 60601-1, IEC 60601-1-2, IEC 60601-2-52

## All Basic American Extended Care Beds Are \_\_\_\_\_

- Certified transport device - passes threshold clearance at ALL heights - ideal for evacuation
- Silver Solutions<sup>8</sup> powder-coating provides antimicrobial surface protection and deters the growth of odor-causing bacteria
- No holes on top surface of tubes that could harbor fluids/bacteria
- Grounded bed with right-angle hospital grade plug for safety
- Adjustable Wallsaver protects the wall with smallest footprint
- Patented Hi-Lo mechanism adds life to the motors
- Sturdy, no-sag, grid sleep surface with cross braces
- Bed provides 360° mobility for cleaning
- Proudly Made in the USA
- Assist devices available\*
- Tool-less assembly



### Adjustable Wallsaver

- Maximum wall protection
- Adjustable clearance lengths
- Standard on all beds

Optional  
Zenith 7200  
Pendant Shown



Slide-W-I-D-E<sup>®</sup>  
Models Only

### Mattress Retainers<sup>9</sup>

- Adjustable for multiple mattress lengths
- Standard on all beds



All Other  
Models

## Zenith & Matrix Series Bed Standard Features \_\_\_\_\_

- Advanced Positioning System (APS) includes one-button comfort/cardiac chair position (optional on Zenith 7200)
- Straight-line Hi-Lo Travel with patent-awarded design technology
- 'Plug & Play' digital electronics enable future additions
- Energy saving electronics
- Sealed ball bearings joints provide continuous, maintenance-free, silent operation
- Intuitive two-pedal locking system - centrally located, easily visible, and accessible at all heights and widths
- Tool-less assembly

## Zenith & Matrix Series Bed Configurations \_\_\_\_\_

- Embedded staff control with safety-inspired, two-stage lockouts enabling caregivers to lock-out features to the resident
- Standard or Advanced Positioning System (APS)
- APS Staff Controls are optional. (APS Staff Control includes one-button comfort/cardiac chair and Trendelenburg/Reverse Trendelenburg positions)
- Hand Pendant to control Trendelenburg/Reverse Trendelenburg (sold separately)
- Accessories and some options are sold separately\*



### Intuitive Two-Pedal Locking System

- Centrally positioned and easily visible from a distance ensures the caregiver that bed is locked at a glance
- Bed rolls at any height increasing caregiver safety
- Readily accessible at ALL bed widths and heights (even lowest)
- Select Option at time of purchase



## Zenith 9200 SERIES

9200 WITH APS SHOWN WITH  
Slide-W-I-D-E® deck option

# Zenith


## Zenith Series Standard Features

- Straight-Line Hi-Lo Travel with Patent-Awarded Design Technology
- Mounted Head and Footboards create a built-in push bar and reduce breakage



## Zenith 7200 SERIES

7200 & 7200 WITH APS

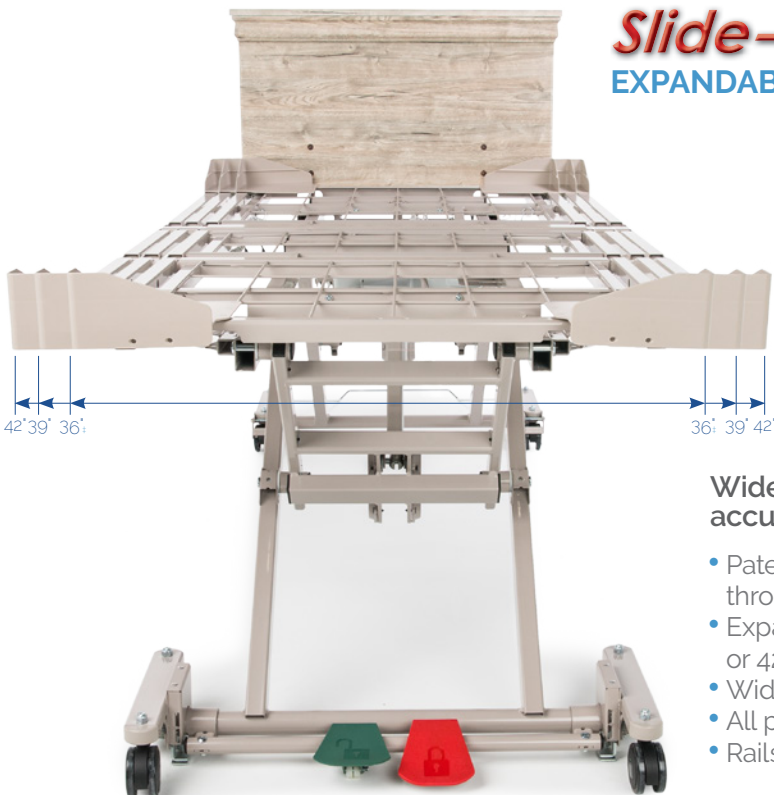

## Slide-W-I-D-E®

EXPANDABLE/RETRACTABLE SLEEP DECK  
OPTIONAL ON ZENITH 7200

Quick-release  
ball-detent pins  
attached to beds  
with lanyards



Wider sleep surfaces may help reduce falls by residents accustomed to sleeping on common consumer products.<sup>7</sup>

- Patented design allows for quick adjustment, ideal for moving through doorways to assist in emergency evacuations
- Expandable/retractable deck quickly and easily adjusts to 36", 39", or 42 inches with two quick-release ball-detent pins on each side
- Wide head and footboards available in 42 inch width
- All pins are attached to the bed with lanyards, eliminating lost parts
- Rails and assist devices are quick-release without tools





**Matrix 6200 APS** SHOWN  
with intuitive pedal-lock option  
Raised Head Position

# Matrix

## Matrix Series Standard Features

- Advanced Positioning System (APS) includes One-Button Comfort/Cardiac Chair Position
- Straight-Line Hi-Lo Travel with Patent-Awarded Design Technology

## Matrix 6200 SERIES 6200 WITH APS

## Advanced Positioning System (APS)

Advanced positioning allows the bed to go into numerous positions

- Positively influencing clinically driven outcomes through multiple positioning options
- Accommodating post-surgical rehab positioning to facilitate resident comfort
- May reduce caregiver strain and potential injuries
- Built-In Hi-Lo Lockout Feature (on APS Hand Controls)



Reverse Trendelenburg




Trendelenburg Position



**Zenith 9200 APS** SHOWN  
with Advanced Positioning System, Slide-W-I-D-E®  
deck option, and intuitive pedal-lock option  
Comfort/Cardiac Chair Position



**Liberty** SHOWN  
with intuitive pedal-lock option  
Supine Position

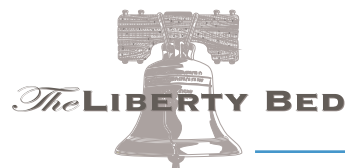

## Liberty Bed Standard Features

- Desired features at affordable price point
- Four Locking casters

## Liberty Bed Available with

- Optional one-step red/green pedal lock to maximize your security (*as shown*)
- Optional 4-inch Total-lock casters increase under-bed clearance

## Overbed Table Options








Rectangle



Raised Panel



Mission

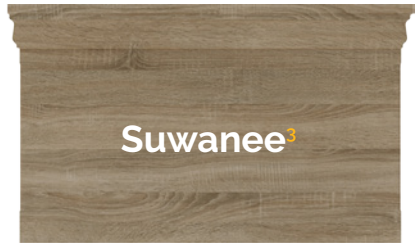

Suwanee<sup>3</sup>



Center Crown




Crescent with Cap<sup>3</sup>



Crescent

3DL material and custom finishes available for all board styles upon request - costs may vary

## Counter-Rotating Assist Device

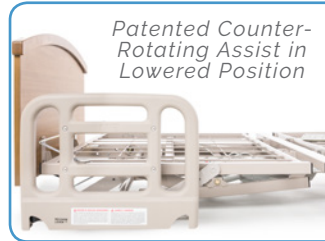


Patented Counter-Rotating Assist in Lowered Position




Raised Position

Provides a sturdy and secure hand hold to assist residents in and out of bed yet easily swivels out of the way when desired

- Meets FDA Entrapment Guidelines
- Patented Counter-Rotating Assist Devices
- Removable without tools
- Moves out of the way for lateral transfers
- Pivots away from egress area
- Provides toe clearance in lowest position to help prevent staff injury
- Hand pendant mounts facing inward or outward and away from hand grab area

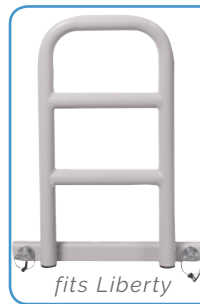


fits Liberty

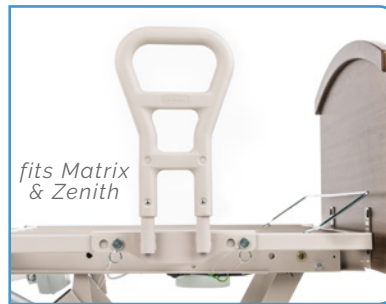

fits Matrix & Zenith

## Pivot Assist Bars

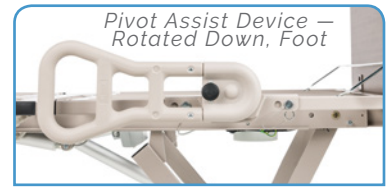

Pivot Assist Device - Rotated Down, Foot



Pivot Assist Device - Up-Right Position



Pivot Assist Device - Rotated Down, Head

## Fixed Assist Devices

Provides a sturdy and secure hand hold to assist residents in and out of bed

- Moves out of the way for lateral transfers
- Tool-free assembly and removal



## Bed Cart

- Ideal for ease of storage
- Allows beds to fit through narrow doorways
- Designed for One-Person Mounting and maneuverability
- Compatible with all Basic American models



## USB Power Supply

- Charges most devices
- Compatible with all Basic American models



## Side Wall Bumper

- Maximizes wall protection
- Compatible with all Zenith and Matrix models

## Advanced Positioning Pendant

- Upgrade from Standard to APS on Zenith 7200
- One-Button Comfort/Cardiac Chair Position
- Hi-Lo Lockout Feature








## Board Mounted IV Socket

- Mounts with 2 screws
- Easy to install
- IV pole sold separately
- Compatible with all Basic American models<sup>2</sup>

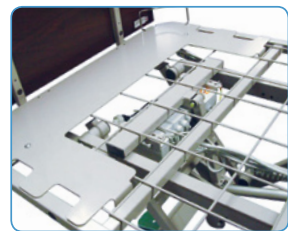

## Surface Length Extension Kit

- Extends standard length bed frame from 84" to 88". Complete with solid pan extension, mattress guard, and installation hardware
- Compatible with all Basic American models<sup>2</sup>
- Easy to install

## UnderBed Light

- Helps reduce fall risks
- Eliminates bright light intrusion
- Multi-level light intensity
- Compatible with all Basic American models



## Embedded Staff Control

- Standard Functionality and Lockouts, or Advanced Positioning (APS)
- Angled-Up Display for Improved Visibility
- Two-Stage lockouts deter unwanted changes
- Mix or Match with Plug & Play design
- Upgrade anytime<sup>2</sup>




## Plastic Pan Deck Kit

- Compatible with all Zenith and Matrix models<sup>2</sup>



## Edema Foot Ratchet Kit

- Elevated foot rest
- Standard on Zenith 9200
- Optional on Matrix 6200 and Zenith 7200

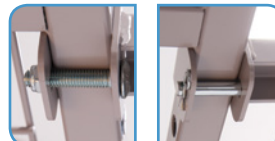

Bolt & Quick Release Pin shown on Wide Kit



## Pin-on Wide Kit

- Compatible with all 36"<sup>+</sup> wide Zenith and Matrix models
- Available in 39" or 42" width kits
- Anytime add-on<sup>2</sup>



## Trapeze and Adapter

- Pivots out of the way when not in use
- Compatible with all Zenith and Matrix models
- Rigid Design

Kensington Maple

Solar Oak

Fonthill Pear

Wild Cherry

Windsor Mahogany

Shaker Cherry

Talas Cherry

Aria

Cocobala

## Bed Specifications:

|   | Zenith 9200                      | Zenith 7200                      | Matrix 6200                | Liberty          |
|---|----------------------------------|----------------------------------|----------------------------|------------------|
| Adjustable Patient Surface Length:  | 80 - 84 inch                     | 80 - 84 inch                     | 80 - 84 inch               | 76 - 80 inch     |
| Fixed Sleep Deck Width:   | 36 <sup>†</sup> inch             | 36 <sup>†</sup> inch             | 36 <sup>†</sup> or 42 inch | 36 inch          |
| Slide-W-I-D-E™ Option Sleep Deck Width:   | 36 <sup>†</sup> or 39 or 42 inch | 36 <sup>†</sup> or 39 or 42 inch | —                          | —                |
| Height Travel Range:  | 7 - 30 inch                      | 7.95 - 30 inch                   | 8.95 - 27 inch             | 9.95 - 24.4 inch |
| Safe Working Load:<br><i>includes resident, support surface, head/footboards, and accessories</i> | 600 lb                           | 550 lb                           | 450 lb                     | 450 lb           |
| Maximum Resident Weight:<br><i>EVENLY DISTRIBUTED - refer to user manual</i>                      | 550 lb                           | 500 lb                           | 400 lb                     | 400 lb           |
| Actuators:  | 4 DC                             | 4 DC                             | 4 DC                       | 3 DC             |

## Options & Features:

|   |         |             |             |                       |
|---|---------|-------------|-------------|-----------------------|
| Chair Position (APS):   | ■       | ■           | ■           | —                     |
| Advanced Positioning Pendant:   | ■       | ■           | ■           | —                     |
| Pendant Attaches to Either Side of Bed:   | ■       | ■           | ■           | —                     |
| Trendelenburg/Reverse Trendelenburg Positioning:<br><i>via APS staff control</i>  | ■       | ■           | ■           | —                     |
| Slide-W-I-D-E® Option Available:<br><i>expands to 36<sup>†</sup>, 39, or 42 inches</i>  | ■       | ■           | —           | —                     |
| Patented Super-Efficient Transmission:<br><i>Improved speed and less strain on the motor</i>                                  | ■       | —           | —           | —                     |
| Hi-Low Actuators are Unloaded in Low Position:  | ■       | ■           | ■           | —                     |
| Embedded Staff Control: <sup>2</sup>  | ■       | ■           | ■           | —                     |
| Upgrade Hand Controls and/or Staff Controls Anytime:<br><i>interchangeable, mix and/or match standard and/or APS controls</i> | ■       | ■           | ■           | —                     |
| Splash Proof Electronics:   | ■       | ■           | ■           | —                     |
| Intuitive Pedal Lock or 4 Locking Casters:  | ■       | ■           | ■           | ■ <sup>2</sup>        |
| Bed Locks and Rolls at Any Height:<br><i>reducing care-giver back strain</i>  | ■       | ■           | ■           | ■                     |
| Silent Operating, Sealed Ball Bearings:<br><i>Maintenance-Free</i>  | ■       | ■           | ■           | —                     |
| Taller Casters:<br><i>approximately 7 inches of clearance</i>   | —       | —           | —           | ■                     |
| Headboard and Footboard Mounting:   | Tubular | Low Profile | Low Profile | Low Profile           |
| Built-in Push Bar and Breakage Reduction:   | ■       | —           | —           | —                     |
| Maximum Clearance:<br><i>overbed tables and lifts</i>   | ■       | ■           | ■           | 5.1 inch <sup>†</sup> |
| Can be Used as Transport Devices:   | ■       | ■           | ■           | —                     |
| FDA Guideline Compliance:<br><i>for Bed Entrapment Risk Reduction<sup>5</sup></i>   | ■       | ■           | ■           | —                     |
| Silver Solutions Powder Coating:<br><i>process provides antimicrobial surface protection<sup>6</sup></i>                      | ■       | ■           | ■           | —                     |
| Made in the USA by Skilled Craftsmen:<br><i>built to withstand trouble-free use for many years</i>                            | ■       | ■           | ■           | —                     |

## Accessories:

|   |   |   |   |   |
|---|---|---|---|---|
| Edema Foot Ratchet: <sup>2</sup><br><i>elevates the foot end of the bed</i> | ■ | ■ | ■ | — |
| Plastic Pan Deck: <sup>2</sup>  | ■ | ■ | ■ | — |
| Backup Battery: <sup>2</sup><br><i>(external or mounted onboard)</i>        | ■ | ■ | ■ | — |
| Extension Kit: <sup>2</sup><br><i>(converts from 84 to 88 inch length)</i>  | ■ | ■ | — | — |
| Pin-on Wide Kit: <sup>2</sup>   | ■ | ■ | ■ | — |

## Warranty:

|                                |          |          |          |          |
|--------------------------------|----------|----------|----------|----------|
| Main Frame & Welds             | 15 years | 15 years | 15 years | 10 years |
| Control Box & Actuator Motors  | 5 years  | 5 years  | 4 years  | 2 years  |
| Hand Control Pendant & Cabling | 3 years  | 3 years  | 3 years  | 1 year   |
| Staff Control & Cabling        | 3 years  | 3 years  | 3 years  | —        |
| Head & Footboards              | 1 year   | 1 year   | 1 year   | 1 year   |

## Product Page Link

|  |  |  |  |
|--|--|--|--|
|  |  |  |  |
|--|--|--|--|

■ Standard   ■ Available   ■ Optional   — Not Applicable/Available

\* options available at additional cost - accessories sold separately   † fixed base   ‡ ± 1 inch tolerance   1 maximum weight capacity - evenly distributed   2 anytime add-on or with purchase   3 boards constructed with 3DL material  
<sup>4</sup> Zenith and Matrix models   <sup>5</sup> Hospital Bed System Dimensional and Assessment Guidance to Reduce Entrapment\* details at www.fda.gov   <sup>6</sup> inhibits the growth of bacterial odors, resists stains caused by mildew and/or mold, preserves the cured coating during the service life of the bed, and inhibits the growth of organisms that can detract from the appearance and coating functionality   <sup>7</sup> Fragala G, Perry B, Fragala M, Examining bed width as a contributor to risk of falls from bed in long term care. Annals of Long-Term Care: Clinical Care and Aging, 2012; 20(6):35-38   <sup>8</sup> GF Health Products, Inc. makes no independent representations or empirical statements about Silver Solutions. This information is reprinted from: TCI Powder Coatings, EPA Guidance for EPA Registration 40 CFR 152.25  
<sup>9</sup> wireframe retainers are standard on all beds except Slide-W-I-D-E® models which use plastic corner retainers