

BARIATRIC/HEAVY DUTY WHEELCHAIRS

TRAVELER® HTC WHEELCHAIR

The Traveler HTC is manufactured with the same commitment to quality E&J is known for. It comes with standard O2 holder with a telescoping IV pole, an attendant push handle bar and rear mounted wheel locks that offer easier maneuvering and control for the caregiver.

Designed for optimal patient transport as well as hospital use, the HTC is built with a heavy duty reinforced frame at all critical load-bearing points. Features include an embossed back and seat upholstery with a double inner liner for durability. Inventory controls are further enhanced by incorporating fixed desk arms, anti-tippers, as well as anti-theft nuts on the front riggings as standard features.

The Traveler HTC provides an all-in-one heavy duty solution with anti-theft control mechanisms along with enhanced patient care features.

Features

- Available in 20", 22" and 24" seat widths, with 18" seat depth (1" longer than standard)
- Carbon steel frame include double cross braces, heavy-duty 5/8" diameter axles
- Embossed black leatherette upholstery with double inner lining for durability with standard chart pocket on back
- Hammertone paint provides a long lasting and durable finish
- Two-position axle and caster fork provide 2" seat height adjustment for a hemi position
- Fixed attendant push bar/handle
- Fixed desk arms and lock nuts on the front rigging ensure inventory controls
- Anti-tippers and one black painted O2/IV telescopic holder for hospital/clinical settings come standard (O2/IV Pole is included as standard after serial number HYT039WC)
- Tool-free adjustable elevating legrests with durable aluminum black footplates
- Rear 24" x 1" black plastic wheels with black handrims and non-marking polyurethane tires
- Maximum Weight Capacity 500 lb -EVENLY DISTRIBUTED
- Meets Medicare code K0007

Warranty:

Five (5) year limited warranty



Rear mounted attendant wheel lock



Push handle



Anti-theft nuts on front riggings

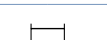



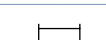























*fixed desk arms (shown) and
elevating legrest (shown)
— other configurations available*





Traveler® HTC Heavy Duty Wheelchair

Item	HCPCS [‡] code	Overall Chair Width		Overall Chair Depth		Overall Chair Height		Folded Chair Width		Seat Width		Seat Depth		Seat Height		Back Height		Wheel Diameter		Caster Diameter		Product [§] Weight w/o front rigging		Weight [†] Capacity EVENLY DISTRIBUTED		Fixed Desk arms	Swingaway Footrest front rigging	Elevating Legrest	Rear Anti-Tippers
																													
		inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	inch	cm	lb	kg†	lb	kg†				
3J010330	K0007	28.5	72.4	34	86.4	36	91.4	—	—	20	50.8	18	45.7	20	50.8	16.5	41.9	24	61	8	20.3	60	27.2	500	226.8		—		
3J010430	K0007	30.5	77.5	34	86.4	36	91.4	—	—	22	55.9	18	45.7	20	50.8	16.5	41.9	24	61	8	20.3	60	27.2	500	226.8		—		
3J010530	K0007	32.5	82.6	34	86.4	36	91.4	—	—	24	61	18	45.7	20	50.8	16.5	41.9	24	61	8	20.3	60	27.2	500	226.8		—		
		 Standard	 Optional	 Available	 Accessory	 Not Apply / Available																							

■ Standard ■ Optional ■ Available ■ Accessory — Not Apply / Available

[‡] designed to meet the requirements of [†] standard rounding method used [†] maximum weight capacity - evenly distributed - refer to user manual
[§] approximate measurement; product weight — as seat width and rigging change, so does product weight ★ options available at additional cost - accessories sold separately

Parts & Accessories: ★ sold separately

	HCPCS [‡]	UOM	Item No.
Anti-Fold Bar (Chromed) Lower Frame	—	each	92523850
Anti-tippers (includes hardware) for HTC	E0971	pair	90763650

	HCPCS [‡]	UOM	Item No.
Positioning Belt (auto-style, push button) [adjustable length: 27 - 47 inches 168.58 - 119.38 cm]	E0978	each	300-505
Attendant wheel lock, rear wheel (includes hardware) for HTC	—	pair	90763640

	HCPCS [‡]	UOM	Item No.
Anti-theft locking nuts for front rigging	—	4 each	90763630
Wheelchair Positioning Tray	—	each	EJTRAY18



GF Health Products, Inc. is an
ISO 13485:2016 / MDSAP Certified Company

Color variations may occur due to monitor or printer settings.

+1 770.368.4700

Information contained herein is subject to change without notice.
The most current and complete product information can be found on our website.

www.grahamfield.com



© 2019, GF Health Products, Inc. All Rights Reserved.
Everest & Jennings is a registered trademark of GF Health Products, Inc.

graham field

GF1200128Rev125