

PATRIOT HOMECARE BED



PATRIOT SEMI-ELECTRIC BED

FEATURES AND ADVANTAGES:

Quality Craftsmanship and Manufacturing

- Designed and engineered by the Extended Care Beds team.
- Sleep surface and motors are **Made in the USA**
- Minimal lead times

SilverSolutions Finish

- Anti-microbial technology that provides protection against bacteria and odor
- Located in the powder-coating providing surface protection
- Easier to clean

Grid Deck Sleep Surface & Frame

- Complete welded steel frame construction makes a much stronger frame with better mattress support
- No spring sag or bottoming out with heavier patients
- The head and foot sections can be raised or lowered simultaneously or independently with one hand control
- Provides better mattress support than traditional spring decks
- Offers maximum pinch point protection between moving deck angles - thus reducing the risk of injury during bed setup
- No connecting links or springs to lose or to connect Grid surface
- Enables one-person delivery and quicker set-up
- Slots in frame ensure quick and proper rail placement
- Assembles easily with just five components

Bed-End Options

- Offered with Fiberboard Decorative Walnut or High-Impact Plastic Bed Ends
- Easy mount bed ends feature locking pins to secure bed ends to sleep deck

Semi-Electric Bed

- The motor unit and bed design provide the ultimate in quiet, smooth operation
- Makes bed height adjustment easy
- Conventional or low bed height for better inventory management!
- Single motor unit is designed for easy installation,
- Requires less maintenance
- During a power failure the motor unit can lower both the head and foot section with just a single standard 9 volt battery

Development and Testing

- Designed to meet the requirements of applicable HCPCS codes
- FDA Hospital Bed System Dimensional Guidance to Reduce Entrapment
- Hospital Bed Safety Workgroup tested
- FDA Dimensional Test Methods for Bed System Zones 1-4 Roll clearance testing
- International Standard IEC60601-2-38 Requirements for the Safety of Electrically Operated Hospital Beds

Options & Accessories*

- The only homecare bed to offer
 - Optional corner mattress retainers keep mattress secure on sleep deck
 - 4 inch extension kit accessory can be added to either the head, foot, or both

 Product Info



Locking pins secure frame to sleep deck

450 POUNDS*
SAFE WORKING LOAD

US0208PL
SEMI-ELECTRIC PATRIOT BED
with High Impact Plastic Bed-Ends
same unit used in cover image

Bed Options.*

Bed-Ends Item Number		HCPCS ¹	Full-Electric Bed	1633 Foam Mattress	Full Rails pair	Half Rails pair	Full Rails pair	Half Rails pair	Quarter Rails pair
Fiberboard Decorative Walnut	High Impact Plastic								
US0208	US0208PL	E0295	■						
US0208-RFPKG	US0208PL-RFPKG	E0260	■	■	■				
US0208-RFPKGHR	US0208PL-RFPKGHR	E0260	■	■		■			
US0208-RPKG	US0208PL-RPKG	E0260	■		■				
US0208-RPKGHR	US0208PL-RPKGHR	E0260	■			■			
US0208-RFPKG	US0208PL-RFPKG	E0260	■	■			■		
US0208-RFPKGHR	US0208PL-RFPKGHR	E0260	■	■				■	
US0208-RFPKGQR	US0208PL-RFPKGQR	E0260	■	■					■
US0208-RPKG	US0208PL-RPKG	E0261	■				■		
US0208-RPKGHR	US0208PL-RPKGHR	E0261	■					■	
US0208-RPKGQR	US0208PL-RPKGQR	E0261	■						■

Warranty

Two year warranty for defects in workmanship and materials of mechanical components, frame, and electronics

During the warranty period, defective items will be repaired or replaced at manufacturer's option. The warranty does not include any labor charges incurred in replacement part(s) installation or any associated freight or shipping charges to the manufacturer

Specifications

	inch	cm [†]	with 4 inch extension kit			
			one end		both ends	
			inch	cm [†]	inch	cm [†]
Overall Width	36.00	91.44	91.00	231.14	95.00	241.30
Sleep Deck Width	36.00	91.44				
Overall Length	87.00	220.98	84.00	213.36	88.00	223.52
Sleep Deck Length	80.00	203.20				
Deck Height Travel Range minimum	15.00	38.10				
Deck Height Travel Range maximum	24.00	60.96				

Electric Height Adjustment	inch	cm [†]
Deck Height Minimum casters in low position	9.50	24.13
Deck Height Maximum casters in low position	20.00	50.80
Deck Height Minimum casters in standard position	13.00	33.02
Deck Height Maximum casters in standard position	23.50	59.69

	degrees [°]
Lift Range - head	60
Lift Range - knee	30

	lb	kg [†]
Maximum Resident Weight ¹ EVENLY DISTRIBUTED	350	158.76
Safe Work Load ²	450	204.12
Product Weight ³ bed and motor only	166	75.30

Cartons Shipped

1. Foot section & motor
2. Head section & Hi/Lo rod
3. Bed-end set
4. Mattress: Compress Rolled for cost effective freight⁴
5. Rail set⁵

Mattress Box Size:
14 x 14 x 40 inches






PATRIOT BED ACCESSORIES

Patriot Bed Rails have been designed to offer a safe and secure bed rail system to provide you with the utmost level of quality to meet your homecare bed needs. Developed using the FDA Hospital Bed System Dimensional Guidance to Reduce Entrapment and tested to the Hospital Bed Safety Workgroup and FDA Dimensional Test Methods for Bed System Zones 1-4 and International Standard IEC60601-2-38 Requirements for the Safety of Electrically Operated Hospital Beds, these new FDA Compliant bed rail systems offer added protection and reduce the risk of entrapment.

Features

- Improved rail design decreases the spacing between the crossbars, thus reducing the risk of entrapment
- New bed rail length closes open areas that can occur during full bed articulation
- Rails mount in fixed rail slots on the bed deck to prevent sliding and moving (shown right)
 - Easy installation saves delivery and setup time and ensures proper rail placement
- Spring-loaded, self-adjusting crossbraces feature safety tabs that prevent accidental disengagement when bumped or during application of side pressure
- Rails are powder-coated with an environmentally-friendly, long-lasting hammertone finish and secure gripping surface
- Heavy-duty steel construction

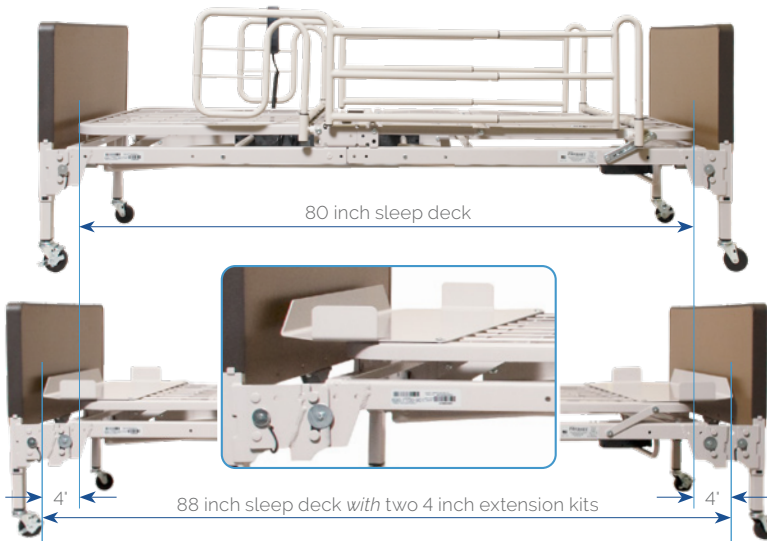

Full Length Rail

The Full Length Rail is a nominal, telescoping full length rail that provides maximum length rail perimeter coverage, intended for use by patients requiring maximum length rail and maximum perimeter protection. Item No: **GF6570B-1**


Bed Extension Kit

Each kit can extend either the head, foot, or both head and foot of the sleep deck 4-inches each up to 88 inches. Includes integrated mattress retainers Item No: **690-0084-000**



Patriot Storage Cart

Bracket set creates a self-contained complete bed transporter that utilizes bed ends to transport the head and foot spring sections, and other bed components. Item No: **600-1003-000**

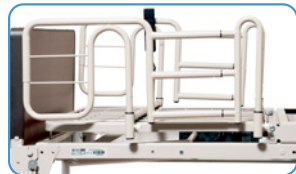


Half No-Gap Head Rail

Provides full head coverage of 31.5 inches and extends to eliminate the gap between the head board and the end of the rail. Item No: **GF6590BH-1**




Quarter Rail

Approximately 18 inches in length. When mounted, this rail can be used as a head assist, head rail, or foot rail to provide minimal rail coverage and gives assistance when entering or exiting the bed. Item No: **GF6580BH-1**

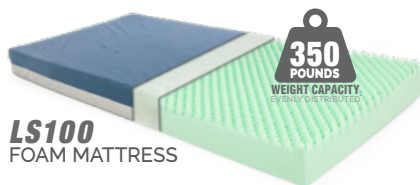

LS100
FOAM MATTRESS

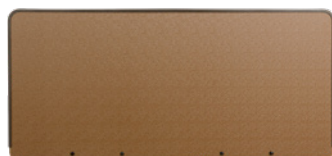

US0208PL
HIGH IMPACT PLASTIC BED-END



US0208
FIBERBOARD DECORATIVE WALNUT BED-END

ALL FUNCTIONS
FULLY-INTEGRATED
INTO CONVENIENT
**HANDHELD
PENDANT**



US0218-RFPKGHR

PATRIOT SEMI-ELECTRIC BED
with Fiberboard Decorative Walnut Bed-Ends and LS100-35 Rolled Foam Mattress and Half Rails

450
POUNDS
SAFE WORKING LOAD

★★★
PATRIOT

★★★ PATRIOT



GF Health Products, Inc. is an
ISO 13485:2016 / MDSAP Certified Company

Images are digital reproductions; perceived texture and coloring may vary on
finished product. Monitor and printer settings will effect appearance of color

•1 770.368.4700

Information contained herein is subject to change without notice. The most
current & complete product information can be found on our website.

www.grahamfield.com



© 2009, GF Health Products, Inc. All Rights Reserved.
Lumex is a registered trademark of GF Health Products, Inc.

published **March 2026**

graham
field

GF09000257RevE26